



Micro-puncher

Convenient Tool for Sampling Film-Like Materials

Simply apply the tip of the tool against a thin or film-like material, and rotate the tool gently, then push out a sample piece.


Variations

id (mm) 0.5 0.75 1.25 2.0 3.0 4.0 5.0



Pushing the button of Micro-puncher, and stainless steel (or a white PTFE) nozzle pushes out a sample piece.

Sampling of papers

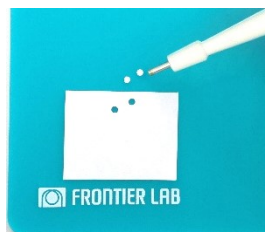

Sampling of leaves



Micro-punchers

Product name	Part number	Specifications	
		Contents	id(mm)
Micro-puncher 050S	FMP-0.50S	puncher (3 pcs), Cutting Mat L (1 sheet)	0.5
Micro-puncher 075S	FMP-0.75S	puncher (3 pcs), Cutting Mat L (1 sheet)	0.75
Micro-puncher 125S	FMP-1.25S	puncher (3 pcs), Cutting Mat L (1 sheet)	1.25
Micro-puncher 200S	FMP-2.00S	puncher (3 pcs), Cutting Mat L (1 sheet)	2.0
Micro-puncher 300S	FMP-3.00S	puncher (3 pcs), Cutting Mat L (1 sheet)	3.0
Micro-puncher 400S	FMP-4.00S	puncher (3 pcs), Cutting Mat L (1 sheet)	4.0
Micro-puncher 500S	FMP-5.00S	puncher (3 pcs), Cutting Mat L (1 sheet)	5.0
Micro-puncher 03 Kit	FMP-03KIT	3 pcs each of 0.75 / 1.25 / 2.0 (mm id) punchers and 2 sheets of Cutting Mat L	
Micro-puncher 04 Kit	FMP-04KIT	3 pcs each of 3.0 / 4.0 / 5.0 (mm id) punchers and 2 sheets of Cutting Mat L	

Cutting Mat

Product name	Part number	Specifications	
		Quantity in package	D x W (mm)
Cutting Mat L	FCM-6577	5	150 x 220

- There is no kit that contains Micro-puncher with 0.5 mm id.
- Sets of 10 Micro-punchers with a same id are available.



FRONTIER LABORATORIES LTD.

www.frontier-lab.com

4-16-20, Saikon, Koriyama, Fukushima, 963-8862 Japan Phone: 81-24-935-5100 Fax: 81-24-935-5102